

**भारतीय मानक**  
**Indian Standard**

**IS 17973 : 2022**  
**ISO 19679 : 2020**

---

---

**प्लास्टिक — समुद्री जल/तलछट  
इंटरफेस में गैर-फ्लोटिंग प्लास्टिक  
सामग्री के एरोबिक बायोडीग्रेडेशन का  
निर्धारण — उद्विकसित कार्बन  
डाइऑक्साइड की विश्लेषण पद्धति  
द्वारा**

**Plastics — Determination of Aerobic  
Biodegradation of Non-Floating  
Plastic Materials in a  
Seawater/Sediment Interface —  
Method by Analysis of Evolved  
Carbon Dioxide**

ICS 83.080.01

© BIS 2022



भारतीय मानक ब्यूरो  
BUREAU OF INDIAN STANDARDS  
मानक भवन, 9 बहादुर शाह ज़फर मार्ग, नई दिल्ली - 110002  
MANAK BHAVAN, 9 BAHADUR SHAH ZAFAR MARG  
NEW DELHI - 110002

[www.bis.gov.in](http://www.bis.gov.in) [www.standardsbis.in](http://www.standardsbis.in)

November 2022

Price Group 8

## Methods of Sampling and Test for Plastics Sectional Committee, PCD 27

### NATIONAL FOREWORD

This Indian Standard which is identical with ISO 19679 : 2020 'Plastics — 'Determination of aerobic biodegradation of non-floating plastic materials in a seawater/sediment interface — Method by analysis of evolved carbon dioxide' issued by the International Organization for Standardization (ISO) was adopted by the Bureau of Indian Standards on the recommendation of the Methods of Sampling and Test for Plastics Sectional Committee and approval of the Petroleum, Coal and Related Products Division Council.

The text of ISO Standard has been approved as suitable for publication as an Indian Standard without deviations. Certain conventions are, however, not identical to those used in Indian Standards. Attention is particularly drawn to the following:

- a) Wherever the words 'International Standard' appear referring to this standard, they should be read as 'Indian Standard'.
- b) Comma (,) has been used as a decimal marker while in Indian Standards, the current practice is to use a point (.) as the decimal marker.

In reporting the result of a test or analysis made in accordance with this standard, if the final value, observed or calculated, is to be rounded off, it shall be done in accordance with IS 2 : 2022 'Rules for rounding off numerical values (*second revision*)'.

# Contents

Page

<b>Introduction</b>	<b>iv</b>
<b>1 Scope</b>	<b>1</b>
<b>2 Normative references</b>	<b>1</b>
<b>3 Terms and definitions</b>	<b>1</b>
<b>4 Principle</b>	<b>2</b>
<b>5 Test environment</b>	<b>2</b>
<b>6 Reagents</b>	<b>2</b>
<b>7 Apparatus</b>	<b>3</b>
<b>8 Procedure</b>	<b>3</b>
8.1 Test material	3
8.2 Reference material	4
8.3 Preparation of the sediment	4
8.4 Test setup	4
8.5 Pre-conditioning phase	4
8.6 Start of the test	5
8.7 Carbon dioxide measurement	5
8.8 End of the test	6
<b>9 Calculation and expression of results</b>	<b>6</b>
9.1 Calculation	6
9.1.1 Amount of CO <sub>2</sub> produced	6
9.1.2 Percentage of biodegradation	8
9.2 Visual inspection	8
9.3 Expression and interpretation of results	8
<b>10 Validity of results</b>	<b>9</b>
<b>11 Test report</b>	<b>9</b>
<b>Annex A (informative) Example of respirometric system based on CO<sub>2</sub> measurement</b>	<b>11</b>
<b>Bibliography</b>	<b>12</b>

**IS 17973 : 2022**  
**ISO 19679 : 2020**

## Introduction

Products made with biodegradable plastics are designed to be recovered by means of organic recycling in composting plants or in anaerobic digesters. The uncontrolled dispersion of biodegradable plastics in natural environments is not desirable. The biodegradability of products cannot be considered as an excuse to spread wastes that should be recovered and recycled. However, test methods to measure rate and level of biodegradation in natural environments (such as soil or the marine environment) are of interest in order to better characterize the behaviour of plastics in these very particular environments. As a matter of fact, some plastics are used in products that are applied in the sea (e.g. fishing gear) and sometimes they can get lost or put willingly in the marine environment. The characterization of biodegradable plastic materials can be enlarged by applying specific test methods that enable the quantitative assessment of biodegradation of plastics exposed to marine sediment and seawater. Plastic products are directly littered or arrive with fresh waters in the pelagic zone (free water). From there, and depending on density, tides, currents, and marine fouling plastics can sink to the sublittoral, and reach the seafloor surface. Many biodegradable plastics have a density higher than 1 and therefore tend to sink. The sediment passes from aerobic to anoxic and finally anaerobic conditions going from the surface (the interface with seawater) into deeper layers, displaying a very steep oxygen gradient.

## Indian Standard

# PLASTICS — DETERMINATION OF AEROBIC BIODEGRADATION OF NON-FLOATING PLASTIC MATERIALS IN A SEAWATER/SEDIMENT INTERFACE — METHOD BY ANALYSIS OF EVOLVED CARBON DIOXIDE

## 1 Scope

This document specifies a test method to determine the degree and rate of aerobic biodegradation of plastic materials when settled on marine sandy sediment at the interface between seawater and the seafloor, by measuring the evolved carbon dioxide (CO<sub>2</sub>). This test method can also be applied to other solid materials.

This test method is a simulation under laboratory conditions of the habitat found in different seawater/sediment-areas in the sea, e.g. in a benthic zone where sunlight reaches the ocean floor (photic zone) that, in marine science, is called sublittoral zone

The determination of biodegradation of plastic materials and other solid materials buried in marine sediment is outside the scope of this document.

**NOTE** Measurement of aerobic biodegradation can also be obtained by monitoring the oxygen consumption, as described in ISO 18830.

The conditions described in this document do not always correspond to the optimum conditions for the maximum degree of biodegradation to occur.

## 2 Normative references

There are no normative references in this document.

## 3 Terms and definitions

For the purposes of this document the following terms and definitions apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- ISO Online browsing platform: available at <https://www.iso.org/obp>
- IEC Electropedia: available at <http://www.electropedia.org/>

### 3.1

#### **theoretical amount of evolved carbon dioxide**

#### **ThCO<sub>2</sub>**

maximum theoretical amount of carbon dioxide evolved after completely oxidising a chemical compound, calculated from the molecular formula or from determination of *total organic carbon (TOC)* (3.2)

Note 1 to entry: It is expressed as mg of carbon dioxide evolved per mg or g of test compound.

### 3.2

#### **total organic carbon**

#### **TOC**

amount of carbon bound in an organic compound

Note 1 to entry: Total organic carbon is expressed as mg of carbon per 100 mg of the compound.

**IS 17973 : 2022**  
**ISO 19679 : 2020**

### 3.3 dissolved organic carbon DOC

part of the organic carbon in water which cannot be removed by specified phase separation methods, for example by centrifugation at 40 000 ms<sup>-2</sup> for 15 min or by membranes with pores of 0,2 µm to 0,45 µm diameter

### 3.4 pre-conditioning phase

pre-incubation of an inoculum under the conditions of the subsequent test in the absence of test material, with the aim to consume potential organic matter present in excess that could disturb biodegradation measurement and to improve the acclimatization of the microorganisms to the test conditions

## 4 Principle

This test method is based on the determination of evolved CO<sub>2</sub> and derives from ISO 14852. The testing medium is based on a solid phase and a liquid phase. The solid phase is a sandy marine sediment laid in the bottom of a closed flask; the liquid phase is a column of natural or artificial sea water, poured on the sediment. The test material is preferably in the form of a film to be laid down on top of the sediment, at the interface between the solid phase and the liquid phase. This is a simulation of an object that has sunk and finally reached the sea floor. The system is contained in a closed flask.

The CO<sub>2</sub> evolved during the microbial degradation is determined by a suitable analytical method. The level of biodegradation is determined by comparing the amount of CO<sub>2</sub> evolved with the theoretical amount (ThCO<sub>2</sub>) and expressed in percentage. The test result is the maximum level of biodegradation, determined from the plateau phase of the biodegradation curve. The principle of a system for measuring evolved CO<sub>2</sub> is given in ISO 14852:2018, Annex A.

The details of interlaboratory testing based on the test method specified in this document are available in Reference [6].

## 5 Test environment

Incubation shall take place in the dark or in diffuse light in an enclosure which is free from vapours inhibitory to microorganisms and which is maintained at a constant temperature, preferably between 15 °C to 25 °C, but not exceeding 28 °C, to an accuracy of ± 2 °C. Any change in temperature shall be justified and clearly indicated in the test report.

NOTE      Temperatures applied in the test can be different from those found in marine environments

## 6 Reagents

**6.1 Distilled or deionized water**, free of toxic substances (copper in particular) and containing less than 2 mg/l of DOC.

### 6.2 Artificial seawater.

Dissolve:

Sodium chloride (NaCl)	22 g
Magnesium chloride hexahydrate (MgCl <sub>2</sub> · 6H <sub>2</sub> O)	9,7 g
Sodium sulfate (Na <sub>2</sub> SO <sub>4</sub> )	3,7 g

Calcium chloride ( $\text{CaCl}_2$ )	1 g
Potassium chloride (KCl)	0,65 g
Sodium hydrogen carbonate ( $\text{NaHCO}_3$ )	0,20 g

in water (6.1) and make up to 1 000 ml.

### 6.3 Natural seawater/sediment.

Take a sample of a sandy sediment and seawater with a shovel beneath the low-water line into a bucket. Transfer the wet sediment together with seawater into sealed containers for transport and fast deliver it to the laboratory. After delivery, conserve the sediment at low temperature (approximately 4 °C) until use. The seawater/sediment sample should be preferably used within 4 weeks after sampling. Record storage time and conditions.

NOTE Seawater and sediment can also be sampled from large, well-running public marine aquaria.

Measure the TOC, pH and nitrogen content of the sediment and of the natural seawater if used instead of artificial seawater. The carbon content of sediment should be in the range of 0,1 % to 2 %.

A preliminary oxidation can be applied to the sediment in order to decrease the organic matter content and the background respiration. Sediment and seawater are fluxed with air and gently stirred (max. 20 r/min to 30 r/min) in a large container for the desired period of time. Include this pre-treatment process in the test report.

## 7 Apparatus

### 7.1 Test flasks.

Biometer flasks of the volume of about 250 ml are appropriate. Reactors with higher volumes can be used, if test conditions are not affected. The vessels shall be located in a constant-temperature room or in a thermostatic apparatus (e.g. water-bath). Stirring can be applied on seawater on condition that it does not disturb the sediment/seawater interface.

NOTE A suitable apparatus is shown in [Figure A.1](#). An example of a stirred apparatus is given in OECD TG 308: 2002, Annex 4<sup>[2]</sup>.

### 7.2 Container for the $\text{CO}_2$ absorber.

A glass beaker to be located in the headspace of the reactor and filled with 10 ml of  $\text{Ba}(\text{OH})_2$  0,0125 mol/l or with 3 ml of KOH 0,5 mol/l.

### 7.3 Analytical balance, shall have a sensitivity of at least 0,1 mg.

### 7.4 pH meter.

## 8 Procedure

### 8.1 Test material

The test material should be in film or sheet form. Cut samples of the test material in the shape of a disk. Disks shall have a smaller diameter than the glass flasks, so that the disks can be easily laid on the bottom of the glass flask.

The sample shall be of known mass and contain sufficient carbon to yield  $\text{CO}_2$  that can be adequately measured by the system used.

**IS 17973 : 2022**  
**ISO 19679 : 2020**

Use a test material concentration of at least 100 mg/l of seawater plus sediment. This mass of the sample should correspond to a TOC of about 60 mg/l. The maximum mass of sample per flask is limited by the oxygen supply to the glass flask. A test material concentration of 150 mg/l to 300 mg/l of seawater plus sediment is recommended.

Calculate the TOC from the chemical formula or determine it by a suitable analytical technique (e.g. elemental analysis or measurement in accordance with ISO 8245) and calculate the  $\text{ThCO}_2$ .

The form and shape of the test material may influence its biodegradation. Similar shapes and thicknesses should preferably be used if different kinds of plastic materials are to be compared.

**NOTE** When the test material in form of film is laid down on the surface of the sediment, it can limit the gas exchange between the water body and the sediment, promoting the formation of anaerobic zones under the test material. In order to reduce this effect, it is possible to perforate the film sample homogeneously over the entire surface.

## 8.2 Reference material

Use ashless cellulose filters as a reference material<sup>1)</sup>. If possible, the TOC, form, and size should be comparable to that of the test material. As a negative control, a non-biodegradable polymer (such as polyethylene) in the same form as the test material shall be used.

## 8.3 Preparation of the sediment

Filter the sediment in a funnel with a coarse filter paper to eliminate excess seawater. Sediment is ready for testing when dripping of sea water is ended. Sediment after filtering is named "wet sediment" hereafter.

## 8.4 Test setup

Provide several flasks, so that the test includes at least the following:

- a) three flasks for the test material (symbol  $F_T$ );
- b) three flasks for the blank (symbol  $F_B$ );
- c) three flasks for reference material (symbol  $F_C$ );
- d) three flasks for negative control (symbol  $F_N$ ).

Two flasks for test material, blank, reference material, and negative control may be used instead of three for screening purposes.

## 8.5 Pre-conditioning phase

In a typical case, use a test flask with a volume of 250 ml. Put 30 g of the wet sediment on the bottom of the flask. Carefully pour 70 ml of natural or artificial seawater. Reactors with higher volumes can be used, if test conditions are not affected. The test should be performed with a water/sediment volume ratio between 3:1 and 5:1 and a sediment layer of about 0,3 cm to 0,5 cm, depending on the granulometry of the sediment.

When using very coarse-grained sediment, the layer can be increased up to 1,5 cm.

Add  $\text{CO}_2$  absorber to the absorber compartments of the test flask in a typical case 3 ml of  $\text{KOH}$  0,5 mol/l or 10 ml of  $\text{Ba(OH)}_2$  0,0125 mol/l. Place the flasks in a constant-temperature environment and allow all vessels to reach the desired temperature. Take the necessary readings and monitor the  $\text{CO}_2$  evolution.

---

1) Laboratory filter paper Whatman n° 42 has been found satisfactory for this purpose and is an example of a suitable product available commercially. This information is given for the convenience of users of this document and does not constitute an endorsement by ISO of this product.



This phase is carried out in order to verify that the endogenous respiration is similar in the different vessels and also to obtain a preliminary oxidation of excess organic matter, in order to start the test with a lower endogenous respiration. The inoculum can be gently stirred in order to accelerate the biodegradation of excess organic matter.

This phase is generally protracted for a week but is possible to extend this time if a high amount of CO<sub>2</sub> evolved is measured.

In case the CO<sub>2</sub> evolution of a vessel is different, reject the diverging vessel or in case of multiple anomalies, start again using new sediment.

## 8.6 Start of the test

Dunk the plastic film sample, cut as described in 8.1, on the sediment of each vessel. Mass of samples (test and reference material) should be about 20 mg each when using a flask with a volume of 250 ml corresponding to an initial test item concentration specified in 8.1. In order to ensure a homogeneous contact between sample and sediment, it is recommended to cover the sample with a suitable cover slip. The cover slips shall also be introduced in blank vessels, for assuring similar conditions.

**NOTE** A suitable cover slip can be made using a common non-biodegradable vinyl-coated fibreglass mosquito net with a fibre diameter of about 280 µm and a 1,8 mm × 1,6 mm mesh.

For an example, see [Annex A](#).

Repeat the procedure for the reference material and the material for the negative control to the respective flasks. Record the mass of the sediment, the sample and the volume of seawater introduced in each vessel.

Nutrients may be supplemented as needed to support microbial diversity and to maintain the capacity to biodegrade the test material. The need and timing of additional nutrients or other appropriate measures may be judged by observation of the temporal course of the biodegradation of the reference substance cellulose. Any addition and the applied method shall be reported in the test report.

## 8.7 Carbon dioxide measurement

**8.7.1** The CO<sub>2</sub> reacts with Ba(OH)<sub>2</sub> and precipitates as BaCO<sub>3</sub>. The amount of CO<sub>2</sub> produced is determined by titrating the remaining Ba(OH)<sub>2</sub> with 0,05 mol/l HCl to a phenolphthalein end-point or by automatic titrator. Because of the static incubation, the BaCO<sub>3</sub> builds up on the surface of the liquid and shall be broken up periodically by shaking the container gently to ensure continued absorption of the evolved CO<sub>2</sub>. This problem can be avoided by using KOH instead of Ba(OH)<sub>2</sub>, which does not form a precipitate.

**NOTE** A discussion on the use of KOH in place of Ba(OH)<sub>2</sub> is reported in Reference [5].

**8.7.2** The containers for the CO<sub>2</sub> absorber shall be removed and titrated before their capacity is exceeded. The period of time will vary with sediments and test materials and increases slowly as the carbon content of the sediment is reduced (a recommended frequency of every 3 to 4 days for the first 2 to 3 weeks and every 1 to 3 weeks, thereafter). At the time of removal of the containers, the reactor should be allowed to sit open so that the air is refreshed before replacing 10 ml of fresh Ba(OH)<sub>2</sub> and resealing the reactor. The reactors should remain open approximately 15 min.

**8.7.3** The CO<sub>2</sub> evolution rate may reach a plateau when all of the accessible carbon has been oxidized. The test may be terminated at this point or earlier, at the discretion of the user. If possible, the residual test material may be extracted from the sediment with an appropriate method and quantified (optional).

**NOTE** The evolved CO<sub>2</sub> can be quantitatively measured also using other suitable methods such as those based on infrared CO<sub>2</sub>-analysers or those based on TOC analysers equipped with an infrared photometer or on gravimetric analysis.

**IS 17973 : 2022**  
**ISO 19679 : 2020**

## 8.8 End of the test

When a constant level of CO<sub>2</sub> evolution is attained (plateau phase reached) and no further biodegradation is expected, the test is considered to be completed. The maximum test period is 24 months. In the case of long test durations, special attention shall be paid to the technical system (e.g. tightness of the test vessels and connections). Any special measures taken, for example, to ensure microbial diversity or to provide sufficient nutrients shall be detailed in the test report. On the last day of the test, measure the pH, acidify all the bottles with 1 ml of concentrated HCl in order to decompose the carbonates and bicarbonates, continue the test for 24 h and finally measure the amount of CO<sub>2</sub> evolved in each of the series of flasks.

## 9 Calculation and expression of results

### 9.1 Calculation


#### 9.1.1 Amount of CO<sub>2</sub> produced

##### 9.1.1.1 Net CO<sub>2</sub> produced

The first step in calculating the amount of CO<sub>2</sub> produced is to correct the test material reactors for endogenous CO<sub>2</sub> production. The control reactor serves as a blank to correct for CO<sub>2</sub> which may be produced through endogenous respiration of the microorganisms. The amount of CO<sub>2</sub> produced by a test material is determined by the difference (in ml of titrant) between the experimental and blank containers. The next step is to convert ml of HCl titrated into mg of CO<sub>2</sub> produced.

##### 9.1.1.2 Ba(OH)<sub>2</sub> used as CO<sub>2</sub> absorber

When CO<sub>2</sub> enters the absorber containers, it reacts in the following manner:

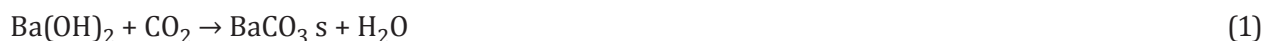

The BaCO<sub>3</sub> formed is insoluble and precipitates. The amount of Ba(OH)<sub>2</sub> remaining in the solution is determined by titration of the 10 ml of the CO<sub>2</sub> absorber with HCl according to the following chemical reaction:



The amount of the remaining Ba(OH)<sub>2</sub> is determined with the [Formula \(3\)](#):

$$R_t = \frac{\text{mol HCl}}{2} \quad (3)$$

where  $R_t$  is the remaining Ba(OH)<sub>2</sub>.

The amount of the reacted Ba(OH)<sub>2</sub> is determined as the difference between the amount originally present in the absorber and the amount remaining after reaction with CO<sub>2</sub>, see [Formula \(4\)](#).

$$R_r = R_0 - R_t \quad (4)$$

where

$R_r$  is the amount of reacted Ba(OH)<sub>2</sub>;

$R_0$  is the amount of Ba(OH)<sub>2</sub> originally present in the absorber;

$R_t$  is the amount of Ba(OH)<sub>2</sub> remaining after the reaction with CO<sub>2</sub>.

This means that the number of mol of CO<sub>2</sub> produced is derived using [Formula \(5\)](#):

$$\text{mol CO}_2 = R_r \quad (5)$$

where  $R_r$  is the amount of reacted Ba(OH)<sub>2</sub>.

### 9.1.1.3 KOH used as CO<sub>2</sub> absorber

The evolved CO<sub>2</sub> reacts with KOH in the following manner:



K<sub>2</sub>CO<sub>3</sub>, the product of [Formula \(6\)](#), is soluble and does not precipitate.

The fresh KOH solution, where no CO<sub>2</sub> has been absorbed, can be titrated with HCl as:



The KOH solutions used as CO<sub>2</sub> absorbers will have both unreacted KOH and K<sub>2</sub>CO<sub>3</sub> as per [Formula \(6\)](#).

During titration, both chemical species will react with HCl, as follows:



The pH shifts in [Formulae \(6\)](#) and [\(7\)](#) are superimposed and cannot be distinguished. Only a single end point in the range of pH between 7 and 8, corresponding to the two reactions, can be identified by using a suitable indicator.

The adsorbed CO<sub>2</sub> can be determined by subtracting from the H<sup>+</sup> equivalents needed to neutralize the original KOH solution and the H<sup>+</sup> equivalents needed to neutralize the reactions represented by [Formulae \(8\)](#) and [\(9\)](#), as shown in [Formula \(10\)](#). In practice:

$$\text{mmol CO}_2 = (V_{\text{HCl}(7)} - V_{\text{HCl}(8+9)}) \times [\text{HCl}] \quad (10)$$

where

$V_{\text{HCl}(7)}$  is the volume of HCl (expressed in ml) consumed in [Formula \(7\)](#);

$V_{\text{HCl}(8+9)}$  is the volume of HCl (expressed in ml) consumed in [Formulae \(8\)](#) and [\(9\)](#);

$[\text{HCl}]$  is the concentration of the HCl solution (0,05 mol/l).

If an end point titrator is available, the mmol of CO<sub>2</sub> can be determined, without using an indicator, with a further reaction. A further addition of HCl makes HCl react with KHCO<sub>3</sub>, produced with [Formula \(9\)](#):



The number of equivalent consumed in [Formula \(11\)](#), and therefore in [Formula \(9\)](#), corresponds to the K<sub>2</sub>CO<sub>3</sub> produced during [Formula \(6\)](#), that in turn corresponds to the absorbed CO<sub>2</sub>.

Consequently, 1 mole of KHCO<sub>3</sub> corresponds to 1 mole of CO<sub>2</sub> reacted in [Formula \(6\)](#), and thus the mmol CO<sub>2</sub> are equivalent to the mmol HCl consumed in [Formula \(11\)](#) end point

**IS 17973 : 2022**  
**ISO 19679 : 2020**

Therefore, the amount of CO<sub>2</sub> can be calculated using [Formula \(12\)](#):

$$mmol\ CO_2 = (V_{HCl(11)}) \times [HCl] \quad (12)$$

where

$V_{HCl(11)}$  is the volume of HCl (expressed in ml) consumed in [Formula \(11\)](#);

$[HCl]$  is the concentration of the HCl solution (0,05 mol/l).

The amount of CO<sub>2</sub> expressed in mg is finally obtained using [Formula \(13\)](#):

$$mg\ CO_2 = mmol\ CO_2 \times 44 \quad (13)$$

where 44 is the molecular weight (g/mol) of CO<sub>2</sub>.

### 9.1.2 Percentage of biodegradation

The percentage of biodegradation is the ration between the evolved CO<sub>2</sub> and theoretical CO<sub>2</sub> (ThCO<sub>2</sub>). The ThCO<sub>2</sub> is shown in [Formula \(14\)](#) and percentage of biodegradation in [Formula \(15\)](#):

$$ThCO_2 = S \times TOC (\%) \times \frac{44}{12} \quad (14)$$

where

$S$  is the amount of specimen (mg)

$TOC (\%)$  is the TOC of the plastic material (or reference material) divided by 100;

44 is the molecular weight (g/mol) of CO<sub>2</sub>;

12 is the molecular weight (g/mol) of C.

Therefore:

$$\%B = \frac{CO_2}{ThCO_2} \times 100 \quad (15)$$

where

$\%B$  is the percentage of biodegradation

CO<sub>2</sub> is the CO<sub>2</sub> produced, expressed in mg;

ThCO<sub>2</sub> is the theoretical amount of evolved CO<sub>2</sub>.

## 9.2 Visual inspection

At the end of the test, check the condition of the samples. If still present, samples can be retrieved for mass determination, other analysis, and photographs.

## 9.3 Expression and interpretation of results

Compile a table of the CO<sub>2</sub> values measured and the percentages of biodegradation for each measurement interval and each test flask. For each vessel, plot an evolved CO<sub>2</sub> cumulative curve and a biodegradation curve in percentage as a function of time.

A curve of mean biodegradation values may be plotted.

The maximum level of biodegradation determined as the mean value of the plateau phase of the biodegradation curve or the highest value, e.g. when the curve decreases or, further on, slowly increases in the plateau phase, characterizes the degree of biodegradation of the test material.

The wettability and the shape of the test material may influence the result obtained, and hence the test procedure may be limited to comparing plastic materials of similar chemical structure.

Information on the toxicity of the test material may be useful in the interpretation of test results showing a low biodegradability.

## 10 Validity of results

The test is considered valid if,

- a) the degree of biodegradation of the reference material ( $F_C$ ) is  $> 60\%$  after 180 days,
- b) the evolved  $\text{CO}_2$  of the blank  $F_B$  at the end of the test does not exceed  $3,5 \text{ mg CO}_2/\text{g wet sediment}$  (see 8.3) after 6 months;

NOTE This value has been determined in an interlaboratory test.

- c) the amount of  $\text{CO}_2$  evolved from the three blank  $F_B$  are within  $20\%$  of the mean at the plateau phase or at the end of the test;
- d) the difference between the percentage biodegradation of the reference material in the different vessels is less than  $20\%$  of the mean at the end of the test;
- e) the percentage of biodegradation of the negative control (flasks  $F_N$ ) is below  $10\%$  at the end of the test.

If these criteria are not fulfilled, repeat the test using another sediment.

## 11 Test report

The test report shall contain at least the following information:

- a) a reference to this document, i.e. ISO 19679:2020;
- b) all information necessary to identify the test and reference materials, including their TOC,  $\text{ThCO}_2$ , chemical;
- c) the main test parameters, including test volume, test medium used, incubation temperature and final pH;
- d) the source and amount of the marine sediment used;
- e) any deviations from the procedure;
- f) any unusual features observed;
- g) the date of the test;
- h) the analytical techniques used, including the principle of the respirometer and the TOC;
- i) all the test results obtained for the test and reference materials (in tabular and graphical form), including the evolved  $\text{CO}_2$ , the percentage biodegradation values;
- j) the duration of the lag phase, biodegradation phase and maximum level of degradation, as well as the total test duration; and, optionally, if run or determined, the negative control  $F_N$ ;

**IS 17973 : 2022**  
**ISO 19679 : 2020**

- k) any other relevant data (e.g. result of the visual final inspection and analysis of final samples, if still retrievable; photos of the final samples);
- l) details of the methods used during the test period in order to support microbial diversity or to avoid nutrient deficiency.


## Annex A (informative)

### Example of respirometric system based on CO<sub>2</sub> measurement

The measurement of the evolved CO<sub>2</sub> can be obtained by trapping the evolved CO<sub>2</sub> in a closed system and then by quantifying it by suitable titration systems. Microorganisms in the vessel consume the oxygen and form CO<sub>2</sub>. This is absorbed by a CO<sub>2</sub> absorber (generally NaOH) that can then be titrated to determine the amount of the absorbed CO<sub>2</sub>.

In a typical case, a 250 ml vessel is used. The sediment occupies about 20 ml, the seawater 70 ml and the headspace 160 ml. The O<sub>2</sub> present in air at 1 atm and 28 °C and a relative humidity of 100 % is about 0,261 mg/ml. This means that the O<sub>2</sub> available at the beginning is 0,261 mg/ml × 160 ml = 41,76 mg (1,305 mmol). The amount dissolved in the seawater can be neglected. This amount of O<sub>2</sub> is sufficient to oxidize to CO<sub>2</sub> and H<sub>2</sub>O an amount of biodegradable organic carbon equal to 15,66 mg and producing 1,305 × 44 = 57,42 mg of CO<sub>2</sub>.

The system is opened in order to refresh the headspace when the O<sub>2</sub> concentration reaches 25 % of the original oxygen concentration.

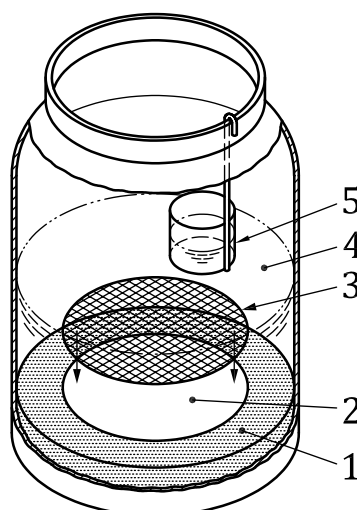

#### Key

- 1 sediment
- 2 sample
- 3 cover slip
- 4 liquid medium
- 5 container for the CO<sub>2</sub> absorber

Figure A.1 — Test flask

**IS 17973 : 2022**  
**ISO 19679 : 2020**

## Bibliography

- [1] ISO 8245, *Water quality — Guidelines for the determination of total organic carbon (TOC) and dissolved organic carbon (DOC)*
- [2] ISO 14852:2018, *Determination of the ultimate aerobic biodegradability of plastic materials in an aqueous medium — Method by analysis of evolved carbon dioxide*
- [3] ISO 18830, *Plastics — Determination of aerobic biodegradation of non-floating plastic materials in a seawater/sandy sediment interface — Method by measuring the oxygen demand in closed respirometer*
- [4] TOSIN M., WEBER M., SIOTTO M., LOTT C., DEGLI-INNOCENTI F., Laboratory test methods to determine the degradation of plastics in marine environmental conditions. *Front. Microbiol.* 2012, **3** p. 225. DOI:10.3389/fmicb.2012.00225
- [5] MODELLI A., CALCAGNO B., SCANDOLA M., Kinetics of Aerobic Polymer Degradation in Soil by Means of the ASTM D 5988 Standard Method. *Journal of Polymer Degradation*, Vol **7**, No. 2, 1999
- [6] ISO/TC 61/SC 5/WG 22 doc. N 17, *Ring test final report seawater sediment interface*, 2015
- [7] OECD TG 308:2002, *Guidelines for the Testing of Chemicals, Test No. 308: Aerobic and Anaerobic Transformation in Aquatic Sediment Systems*



This Page has been Intentionally left blank

## Bureau of Indian Standards

BIS is a statutory institution established under the *Bureau of Indian Standards Act, 2016* to promote harmonious development of the activities of standardization, marking and quality certification of goods and attending to connected matters in the country.

### Copyright

BIS has the copyright of all its publications. No part of these publications may be reproduced in any form without the prior permission in writing of BIS. This does not preclude the free use, in the course of implementing the standard, of necessary details, such as symbols and sizes, type or grade designations. Enquiries relating to copyright be addressed to the Head (Publication & Sales), BIS.

### Review of Indian Standards

Amendments are issued to standards as the need arises on the basis of comments. Standards are also reviewed periodically; a standard along with amendments is reaffirmed when such review indicates that no changes are needed; if the review indicates that changes are needed, it is taken up for revision. Users of Indian Standards should ascertain that they are in possession of the latest amendments or edition by referring to the website- [www.bis.gov.in](http://www.bis.gov.in) or [www.standardsbis.in](http://www.standardsbis.in).

This Indian Standard has been developed from Doc No.: PCD 27 (18151).

### Amendments Issued Since Publication

Amend No.	Date of Issue	Text Affected

## BUREAU OF INDIAN STANDARDS

### Headquarters:

Manak Bhavan, 9 Bahadur Shah Zafar Marg, New Delhi 110002

Telephones: 2323 0131, 2323 3375, 2323 9402

Website: [www.bis.gov.in](http://www.bis.gov.in)

### Regional Offices:

	Telephones
Central : 601/A, Konnectus Tower -1, 6 <sup>th</sup> Floor, DMRC Building, Bhavbhuti Marg, New Delhi 110002	{ 2323 7617
Eastern : 8 <sup>th</sup> Floor, Plot No 7/7 & 7/8, CP Block, Sector V, Salt Lake, Kolkata, West Bengal 700091	{ 2367 0012 2320 9474
Northern : Plot No. 4-A, Sector 27-B, Madhya Marg, Chandigarh 160019	{ 265 9930
Southern : C.I.T. Campus, IV Cross Road, Taramani, Chennai 600113	{ 2254 1442 2254 1216
Western : Plot No. E-9, Road No.-8, MIDC, Andheri (East), Mumbai 400093	{ 2821 8093

**Branches :** AHMEDABAD. BENGALURU. BHOPAL. BHUBANESHWAR. CHANDIGARH. CHENNAI. COIMBATORE. DEHRADUN. DELHI. FARIDABAD. GHAZIABAD. GUWAHATI. HIMACHAL PRADESH. HUBLI. HYDERABAD. JAIPUR. JAMMU & KASHMIR. JAMSHEDPUR. KOCHI. KOLKATA. LUCKNOW. MADURAI. MUMBAI. NAGPUR. NOIDA. PANIPAT. PATNA. PUNE. RAIPUR. RAJKOT. SURAT. VISAKHAPATNAM.

Published by BIS, New Delhi